



Operation Instructions for **CALIDOSE® DISPENSER** for use with Lantheus Xenon Xe-133 Gas

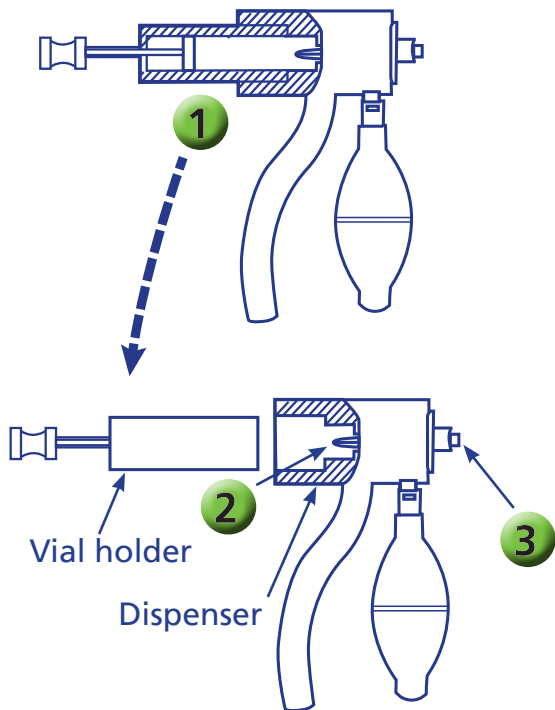

LOADING

1. Separate shielded vial holder from dispenser
2. Check Huber point needles inside the dispenser to ensure that they are not blocked (if necessary, clean by pushing a fine wire through needles)
3. Attach a hypodermic needle (or other connector) securely to Luer Lock fitting on front end of dispenser.
4. Remove cap of the Xenon Xe-133 lead shipping tube, being careful to point opened tube away from body.
5. Place the open end of the shielded vial holder tightly against the open end of the shipping tube.
6. Simultaneously invert both to allow one vial to slip from shipping tube into shielded vial holder.
7. Carefully separate the pieces shown in figure, being certain to point the open ends away from the body. There should be one vial in the shielded vial holder, septum end out. Recap the lead shipping tube.

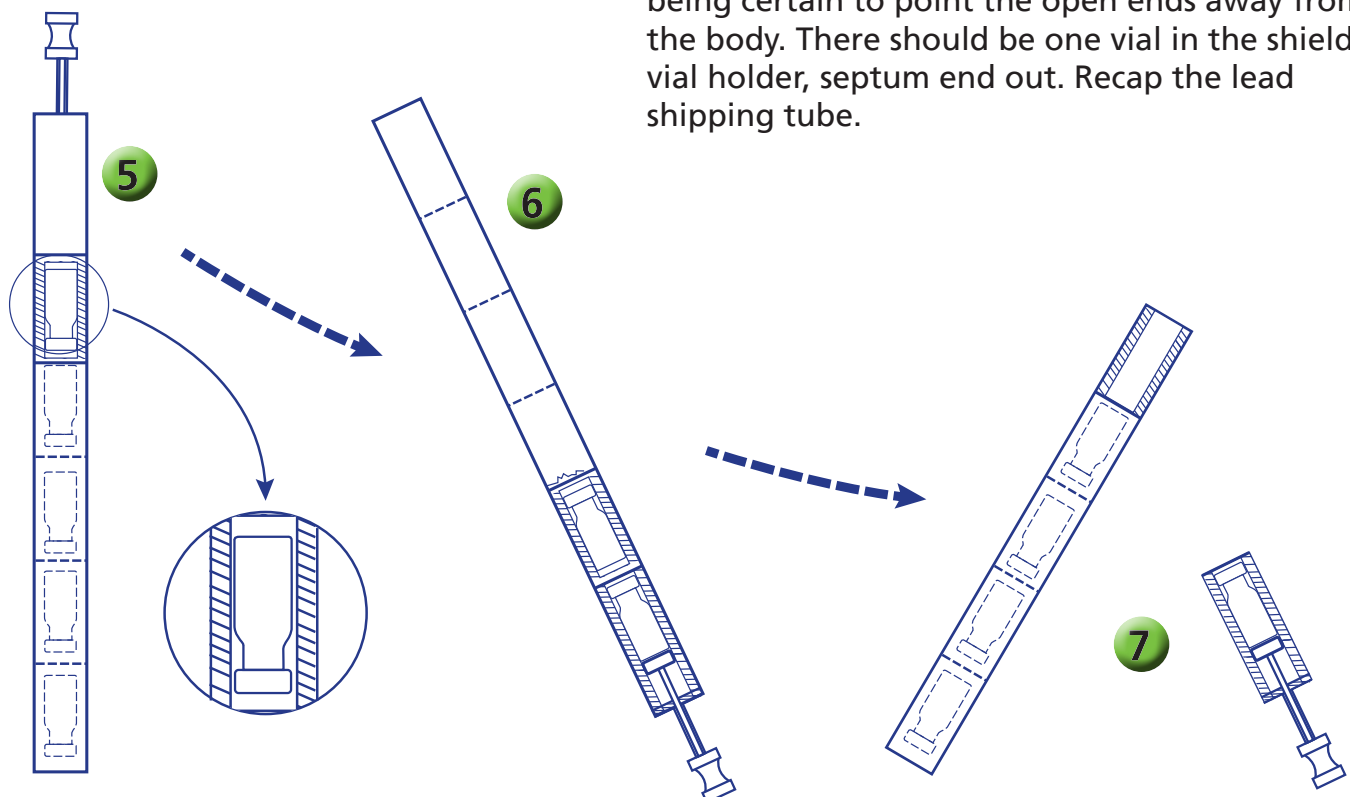

Operation Instructions for **CALIDOSE® DISPENSER** for use with Lantheus Xenon Xe-133 Gas



LOADING (continued)

8. Holding the shielded vial holder upright (and pointed away from the body), insert it into the dispenser until seated. **DO NOT PUSH PLUNGER AT THIS TIME.**

The loaded CALIDOSE® Dispenser is now ready for use and can be put aside until needed. Note that the appropriate radiation protection precautions should be taken with this unit.



USING

1. Affix the CALIDOSE® dispenser to a spirometer or related breathing apparatus.
2. Puncture septum of loaded vial by pushing plunger into the dispenser.
3. Immediately squeeze the rubber bulb, and then release
4. Detach CALIDOSE® assembly from breathing apparatus.

STORING

Remove vial holder from dispenser. The used vial should no longer contain enough residual Xenon Xe-133 to cause any additional radiation exposure to the operator. Once confirmed that the vial has been emptied of its contents, it may be removed by hand for disposal as low level radioactive waste. Replace vial holder in dispenser for easy storage.

The CALIDOSE® Dispenser is for use **ONLY** with Lantheus Xenon Xe-133 vials.

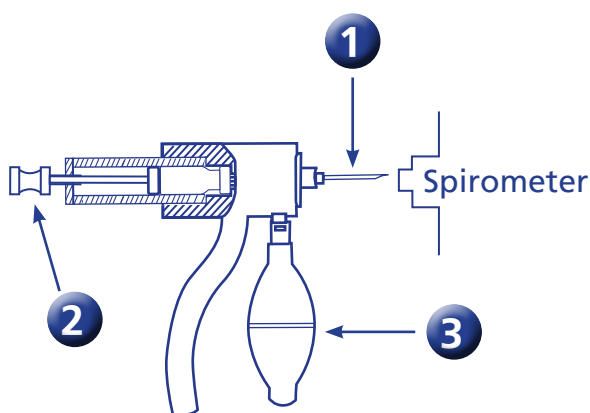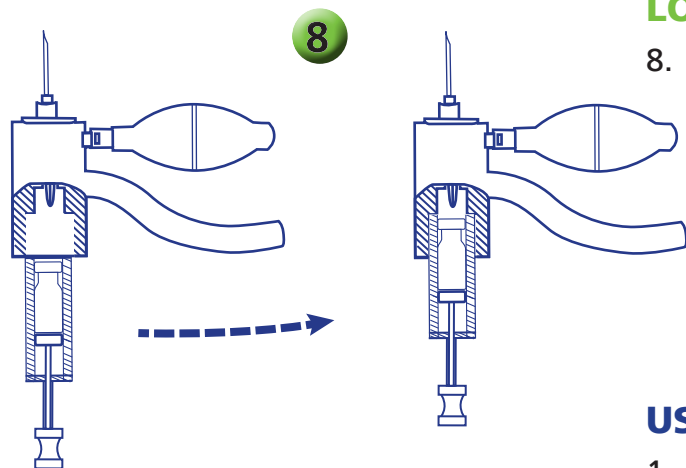

Specify **Lantheus Xenon Xe-133**
when ordering from your local Radiopharmacy

To order from Lantheus, contact your
Lantheus Sales Representative
or call Customer Service: 1-800-299-3431
www.lantheus.com



LANTHEUS™